

Ephedrine and Pseudoephedrine on Lux 3µm AMP 150x4.6mm

Column: Lux® 3 µm AMP, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4,6 mm ID

Order No: 00F-4751-E0

Elution Type: Isocratic

Eluent A: 5mM Ammonium Bicarbonate [pH11]

Eluent B: Methanol

Eluent C: Acetonitrile

Eluent D: Isopropanol

Gradient	Step No.	Time (min)	Pct A	Pct B	Pct C	Pct D
----------	----------	------------	-------	-------	-------	-------

Profile:	1	0	65	20	10	5
----------	---	---	----	----	----	---

Flow Rate: 0,5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)

Analyst Note: Resolution = 1.92, 1.92, 1.72


Selectivity = 1.16, 1.15, 1.13



Products used in this application:



23670

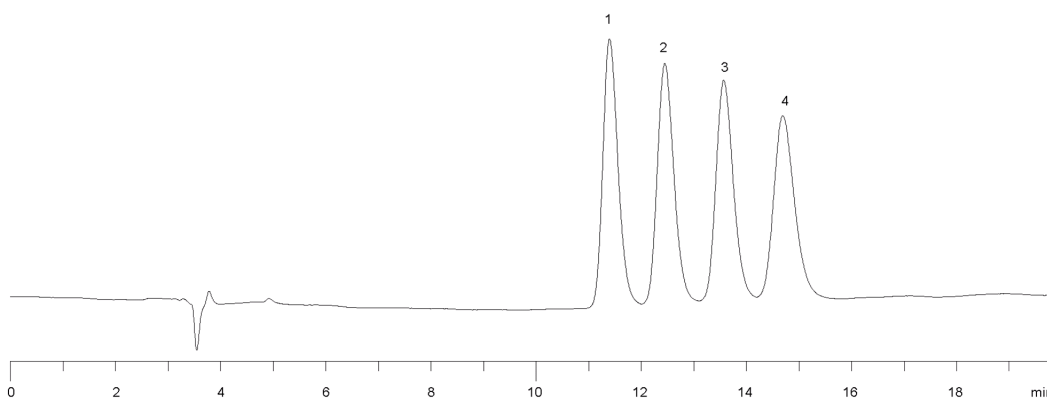

ANALYTES:

1 1S,2R (+)-Ephedrine HCl

Retention Time: 11,39 min

2 1R,2S (-)-Ephedrine HCl

Retention Time: 12,44 min

3 R,R (-)-Pseudoephedrine

Retention Time: 13,58 min

4 S,S (+)-Pseudoephedrine

Retention Time: 14,69 min

