

BTEX in Under Six Minutes on ZB-WAXColumn: Zebtron[™] ZB-WAX, GC Cap. Column 30 m x 0.32 mm x 0.50 μ m, Ea

Phase: Polyethylene Glycol (PEG)

Dimensions: 30 meters x 0,32 mm x 0,5 μ m

Order No: 7HM-G007-17

Oven Profile: 60°C to 75°C at 15°C/min to 90°C at 3°C/min for 3 min


Carrier Gas: Constant Flow Helium, 2 mL/min

Injection: split 20:1 0,2 μ L @ 250°C

Detection: Flame Ionization (FID) (300°C)



Products used in this application:

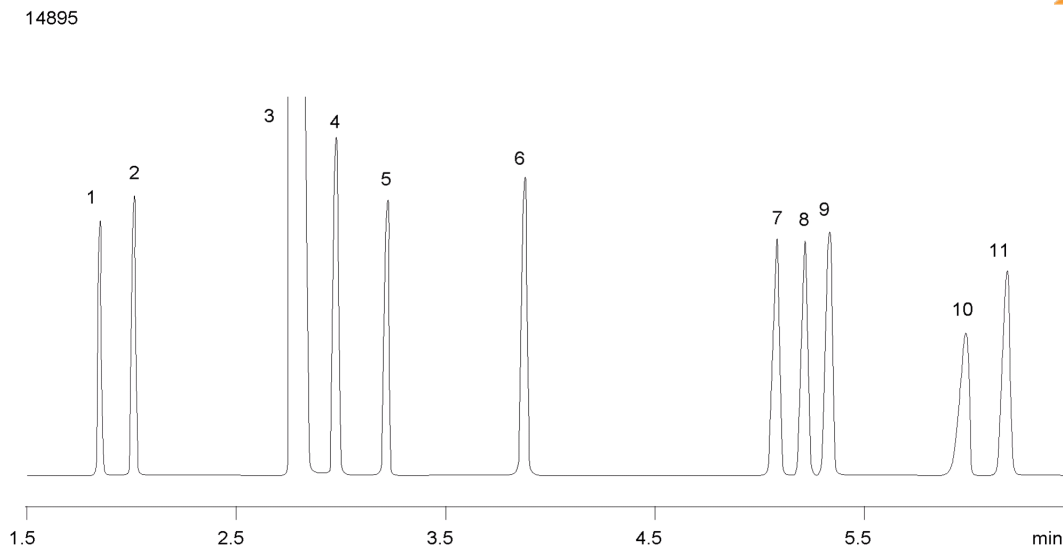

BTEX in Under Six Minutes on ZB-WAX

ANALYTES:

- 1 Pentane
- 2 Heptane
- 3 Methylene chloride (solvent)
- 4 Benzene
- 5 Decane
- 6 Toluene
- 7 Ethylbenzene
- 8 p-Xylene
- 9 m-Xylene
- 10 Dodecane
- 11 o-Xylene

